

Metals Analyses of Portland and Blended Cement Via Microwave Digester and Inductively Coupled Plasma Technologies

Abstract

The method involves a sample preparation technique through microwave digestion to extract eight metals (Si, Al, Mg, Ca, Fe, Na, K, and Ti) in Portland (SRM 1888b) and Blended (SRM 1881a) cement traceable to NIST followed by spectrometric detection and quantitation via the simultaneous system inductively coupled plasma-optical emission spectroscopy (ICP-OES). A 0.2-gram weight of the materials was prepared for closed-system digestion at 1200W power. After digestion, solutions were diluted using 50-mL dedicated polypropylene laboratory wares (volumetric flasks) and were analyzed. Accordingly, the dilution and calibration patterns were determined. The analyses were performed in five replicates. The detections of Si, Al, Mg, Ca, Fe, K, and Ti following radial mode and Na via axial torch orientation gave satisfactory results. Reproducibility studies were conducted, and ICP-OES configurations were established for each element. The method was developed based on its figures of merit. The procedure was validated with exceptional performance except that it was not robust at a higher sample mass. Operational maintenance procedures were established when many samples were run. Measurement uncertainty estimation is recommended to identify the sources and effects of errors on the measurand.

Keywords: Portland and Blended cement; ICP-OES; radial mode; axial; reproducibility studies; figures of merit

I. Introduction

Cement is made by heating a mixture of materials such as chalk or marl, shells, limestone, clay, shale, blast furnace slag, slate, iron ore, and silica sand. The resulting product is a fine powder mainly composed of oxides of silicon, calcium, magnesium, iron, and aluminum, along with small amounts of other elements like potassium and sodium. Cement is valuable because it acts as a binding agent. When mixed with rocks or gravel and combined with water, the mixture hardens into concrete. During hydration, compounds like SiO_2 , CaO , Fe_2O_3 , and Al_2O_3 react to form aluminoferrite hydrates and calcium aluminosilicates, which are responsible for cement's hardening. Two reference materials were used in this validation study, one the Portland (SRM 1888b) and one Blended (SRM 1881a) which is a Portland cement blended with slag and fly ash which when combined, improves the workability of plastic concrete, and the strength and durability of hardened concrete. With such a combination, cost savings and strength increase are benefited.

The performance of cement depends on its mineral composition and its exact chemical composition. For instance, elevated levels of MgO can lead to weakness in the hardened cement

[1]. The characteristics of the raw materials used in cement production can lead to significant variations in the finished product's chemical composition. As a result, chemical analysis of cement is essential for evaluating its quality and determining the specific applications for which a particular cement type is most suitable. Cement testing plays a vital role in ensuring product quality, assessing performance, maintaining quality control, meeting regulatory standards, and promoting cost-effective construction practices. Proper testing helps ensure that concrete structures are safe, durable, and compliant with required quality standards.

Chemical analysis is routinely applied in cement manufacturing for quality control, evaluation of quarry materials, and assessment of intermediate products. Professionals in construction and manufacturing may also analyze cement to confirm its suitability for intended applications and to avoid costly repairs or failures in the future.

Elemental analysis of cement can be used to estimate the proportion of each element present in a sample. Although mineral composition is often more critical in determining cement properties and quality, it is generally more difficult to measure directly.

Compared with mineral analysis, elemental analysis is easier to perform, and the resulting elemental ratios can be converted into weight fractions of the corresponding oxides in the sample. From these oxide ratios, the mineral composition of cement can then be estimated.

A wide range of techniques is available for elemental analysis, including atomic spectroscopy, gravimetric methods, X-ray photoelectron spectroscopy, and inductively coupled plasma (ICP) systems [2,3]. Inductively coupled plasma-optical emission spectroscopy (ICP-OES) instruments can be configured as radial or axial views. In the radial configuration, the plasma is viewed from the side, while in the axial configuration, the plasma is viewed end-on (along the length of the plasma). An axially viewed plasma is used for analysis, typically for samples that require high sensitivity. The limitations are low matrix tolerance and dynamic range resulting from physical effects and the configuration of the plasma torch and interface. For high matrix tolerance, a radially viewed plasma is used in this configuration, and the sensitivity is reduced, but matrix tolerance and dynamic range can be improved. And so, in the case of cement samples, these two modes were tested to assess the sensitivity of elemental components present. However, since extraction or sample preparation is considered the soul of any experiment, the performance capability during microwave digestion was substantial.

All the developed analytical methods need to be validated to make sure that each measurement of the content of the analyte in the sample in the routine analysis is close to the true value [4]. The method was developed based on its reliability characteristics. These include precision, accuracy, sensitivity, linear dynamic range, detection limit, selectivity, and ruggedness and robustness.

II. Experimental

A. Standard and Sample Preparations

Materials and reagents used include pure metal standards (Si, Al, Mg, Ca, Fe, Na, K, Ti) purchased from Merck Millipore. Deionized water, concentrated HNO₃, concentrated HCl, concentrated HF, and concentrated H₂SO₄ were prepared for the digestion of samples. The certified reference materials (CRMs) were SRM 1888b Portland Cement and Blended (SRM 1881a) traceable to NIST. The Anton Paar Multiwave 5000 microwave digester and Shimadzu ICPE-9800 Inductively Coupled Plasma systems were exploited by running replicated samples. Preparation for digestion at 1200W power of both Portland and Blended SRM was made by weighing 200.0 ± 0.1 mg of previously homogenized cement directly into the liners (TFM vessels) of the microwave digester. Actual weights were recorded for recovery calculations. Six milliliters concentrated HNO₃, two milliliters concentrated HCl and 0.5 milliliters HF were added into each liner. Samples were allowed to stand for several minutes with occasional swirling to release some gases until the reaction subsided. Each vessel was sealed with a vent closed, following the recommended program for digestion.

Table 1. Microwave Program

Step	Time	Temperature	Microwave power
1	10 minutes	220°C	Up to 1200 Watt
2	20 minutes	220°C	Up to 1200 Watt

The below figure shows the temperature profile of digestion. As can be seen, a come-up time of 10 minutes was set for the temperature to reach 220°C while heating. A ramping rate of about 19.5°C/minute takes place before digestion at a constant temperature for another 20 minutes. These conditions were achieved using a 1200-watt power.

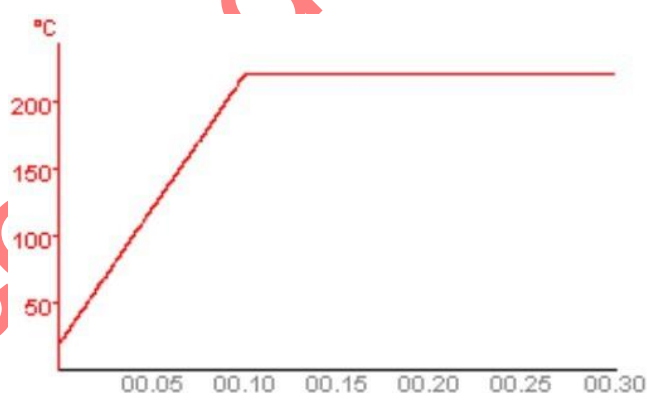


Fig. 1. Temperature profile during sample digestion.

When running more than eight vessels, 1200W power was used and balance inside the cavity was observed. After digesting and cooling, the venting nozzle was opened to release some pressure vent. Digestates were transferred quantitatively into plastic polypropylene 50mL volumetric flasks and made to volume with deionized water (solution A). A dilution pattern was followed accordingly (Table 2).

Table 2. Dilution Pattern (Portland and Blended Cement)

Elements	steps	Dilution rate
Ca, Mg, Al	as is (solution A)	50
Si, Fe, Na, K, Ti	take 5 mL from solution A and dilute to 50 mL	500

Preparation of working standards involved the use of stock solutions of pure metals. Following the calibration pattern below, the volume of stock solutions necessary to prepare a multi-element mixed standards solution was calculated. A calibrated digital pipettor was used for this preparation to guarantee precise and accurate results as well as the linearity of the curves.

Table 3. Calibration Range

Elements	Concentration, ppm		
	I	II	III
Ca	50.0	100.0	200.0
Mg	2.5	5.0	10.0
Si	10.0	20.0	40.0
Al	2.5	5.0	10.0
K	25.0	50.0	100.0
Fe	25.0	50.0	100.0
Na	5.0	10.0	20.0
Ti	5.0	10.0	20.0

B. Quantification with ICP-OES

To operate the ICP, a sequence of steps was conscientiously followed, noting that what comes first during the turning on comes last during the shutting off. First, the breaker was switched on followed by the cooling system ensuring the pump was working. Then, the ICP system was turned on providing a waiting time of at least 4 hours to achieve enough vacuum. The software was opened, and the "Analysis" button was clicked. After that, the argon supply was turned on to ensure the flow rate was properly set. Then, the appropriate method developed was downloaded after which plasma was turned on until all necessary initialization settings were made. Next, the wavelength calibration and registration of all samples and standard information were completed just before the analysis started. After the run, all the necessary data and method information were saved, and data processing and evaluation came next. For data processing, postrun (not real-time) was opened. Selected calibration data files were opened and evaluated to generate calibration curves. Once ready, the processing and integration of samples were completed.

A simple method for microwave-assisted preparation of cement reference samples is shown below.

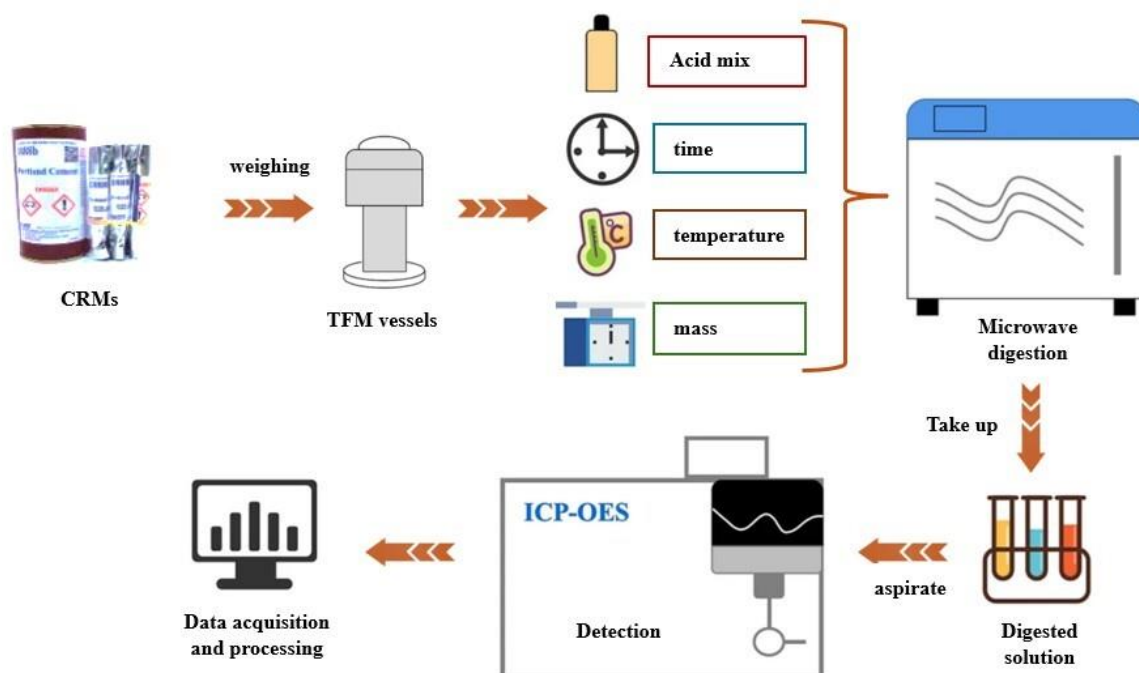


Fig. 2. The microwave-assisted preparation and ICP-OES digestion of cement samples.

C. Wavelength Selection

Each element has a unique spectrum because it possesses a distinct number of protons and a specific arrangement of electrons. These spectral differences arise from variations in the energy absorbed or emitted by atoms when their electrons transition between different energy levels [5]. ICP-OES is an emission spectrophotometric technique in which excited electrons release energy at specific wavelengths as they return to their ground state. A key feature of this process is that each element emits radiation at characteristic wavelengths unique to its chemical identity [6]. However, when the excited electron returns to the ground state, it can take multiple pathways, each corresponding to the emission of a photon of a specific energy. Thus, one absorption can lead to numerous emissions, therefore, more emission lines are observed in the atomic spectra of most elements. Although each element can emit energy at multiple wavelengths, ICP-OES typically focuses on one or only a few selected wavelengths for each element [7]. This approach enhances selectivity in detection. The intensity of the emitted radiation at the chosen wavelength is directly proportional to the concentration of the element in the sample. This can be seen in the generated calibration curves where the linearity and the dynamic range can be measured. ICP sensitivity is in the low parts-per-billion range [8]. The system uses a large-scale CCD (charged coupled detector) with an array of pixels each specific for a particular wavelength, thus enhancing the capture of photons emitted after absorption or during the relaxation stage. That means that the entire chip, which is photoreactive, enables ICP to detect the various emission lines with great accuracy and sensitivity.

D. Reliability Characteristics

A comparison of the amount found versus those of the target declaration on the certificates provides the recovery of the method, which is an estimate of the accuracy of the technique [9]. After the detection and quantitation of the elements of interest, gravimetric conversion factors were used to convert each element into its stoichiometric oxides as components of cement (Table 4). Since the validation was carried out using the CRMs, accuracy was determined by comparing the measured concentration with the certified values and expressed as percentage recovery R (%). This represents the "trueness" or the closeness of the analytical result to the "true" value. The smaller the recovery percentage, the larger the bias of the method thus lower trueness. This measures how the technique can be improved and optimized, from sampling, and preparation to detection. The acceptable recovery percentage depends on the analyte concentration. The acceptability criteria for recovery varies with concentration. The lower the concentration, the lower the probability of getting 100% of the element of interest and the greater the chances of physical and chemical losses (Table 5).

Table 4. Conversion Factors of Each Element to Its Stoichiometric Oxide

Element	Stoichiometric oxides	Ratio	Conversion factor
Fe	Fe ₂ O ₃	$\frac{Fe2O3}{Fe}$	1.4297
K	K ₂ O	$\frac{K2O}{K}$	1.2046
Ti	TiO ₂	$\frac{TiO2}{Ti}$	1.6681
Na	Na ₂ O	$\frac{Na2O}{Na}$	1.3480

Element	Stoichiometric oxides	Ratio	Conversion factor
Ca	CaO	$\frac{CaO}{Ca}$	1.3992
Mg	MgO	$\frac{MgO}{Mg}$	1.6582
Si	SiO ₂	$\frac{SiO2}{Si}$	2.1392
Al	Al ₂ O ₃	$\frac{Al2O3}{Al}$	1.8895

In terms of precision, only reproducibility precision was evaluated. Repeatability precision is the agreement of results from at least seven replicates under common conditions from extraction up to aspiration in the ICP and to assess the routine performance of the method. However, repeatability precision was not evaluated. The normalized intensities of the specific element were one of the bases for choosing the optimum detection wavelength. While the standard deviation of analytical data may remain relatively constant over a narrow range, it generally depends on the magnitude of the values. That is, larger measurements tend to produce larger standard deviations. Therefore, when comparing variability among such data (e.g., assessing precision), it is often more practical to use the relative standard deviation (RSD) rather than the absolute standard deviation. The RSD can be expressed as a fraction, but it is more commonly reported as a percentage, known as the coefficient of variation (CV or %RSD). The usual range for %RSD is 0 – 2.0% for analytical testing. Conversely, the reproducibility precision evaluates the closeness of results from several determinations done under varied conditions [10]. The data set will explain the method's performance across a certain period. This could be based on consensus values among laboratories performing the tests. The proficiency test using the CRMs can provide consensus values and thus be used to evaluate reproducibility precision. Reproducibility precision for each element was determined after getting the mean concentration values from data generated at different conditions. The concentration of the analyte, C, was expressed as a mass fraction, where the numerator and

denominator are in the same mass units, after which the HorRat_R from the RSD_R was evaluated [11]. The predicted reproducibility relative standard deviation calculated from the Horwitz formula:

$$\text{PRSD}_R = 2C^{-0.15} \quad [\text{Eq. 1}]$$

Acceptable reproducibility values typically range from one-half to twice the calculated value. Alternatively, a ratio can be calculated of the found value for RSD_R to that calculated from the formula designated as HorRat_R . Acceptable values for this ratio are typically 0.5 to 2.

$$\text{HorRat}_R = \text{RSD}_R (\text{found, \%}) / \text{PRSD}_R (\text{calculated, \%}) \quad [\text{Eq. 2}]$$

This equation transformed the RSD found to a fraction of the RSD expected from applying the exponential equation to the concentration. The precision found was expressed as a fraction of the precision calculated from the experimental concentration estimates. It equals 1 for exact correspondence. The precision is better than expected if the ratio is less than 1 and poorer if greater than 1 [10]. Table 5 shows the recommended accuracy and reproducibility precision limits as a function of the concentration [12].

Table 5. Acceptable Recovery and Predicted Relative Standard Deviations Criteria

Concentration	Recovery (%)	Reproducibility RSD_R , (%)
100%	98-101	2
10%	95-102	3
1%	92-105	4
0.1%	90-108	6
0.01%	85-110	8
1 $\mu\text{g/g}$ (1 ppm)	80-115	11
1 $\mu\text{g/g}$ (1 ppm)	75-120	16
10 ng/g(10 ppb)	70-125	32

AOAC Guidelines for Single-Laboratory Validation of Chemical Methods for Dietary Supplements and Botanicals (2012) Official Methods of Analysis, 19th Ed., Appendix K, AOAC INTERNATIONAL, Gaithersburg, MD

The linear dynamic range refers to the range of concentrations an instrument can read, from the minimum to the maximum detectable (Figure 3). The range gives satisfactory reliability using the Coefficient of Determination (R^2) values as explained through the next criterion, linearity. Reliability can be judged in terms of the linearity of the curve. Calibration curves must show at least a 0.999 correlation coefficient R value [13]. In other words, R^2 shows how well the data fit the regression model (the goodness of fit).

The limit of detection (LOD) was estimated using the standard deviation of a calibration curve by using the residual standard deviation (s_r), which measures the scatter of the data points around the regression line [14]. The STEYX function in Microsoft Excel directly calculates the residual standard deviation, also known as the standard error of the estimate. It requires the input of the range of the y values and the range of the x values. For LOD, the formula is,

$$\text{LOD} = 3.3 \times s_r / S$$

where, s_r – is the residual standard deviation and

S – is the slope of the calibration curve

This method relies on a linear regression of your calibration data to find the slope and residual standard deviation. To calculate the Limit of Quantification (LOQ), the formula below was used.

$$LOQ = 10 \times s_r/S$$

Figures of merit are described in Figure 3.

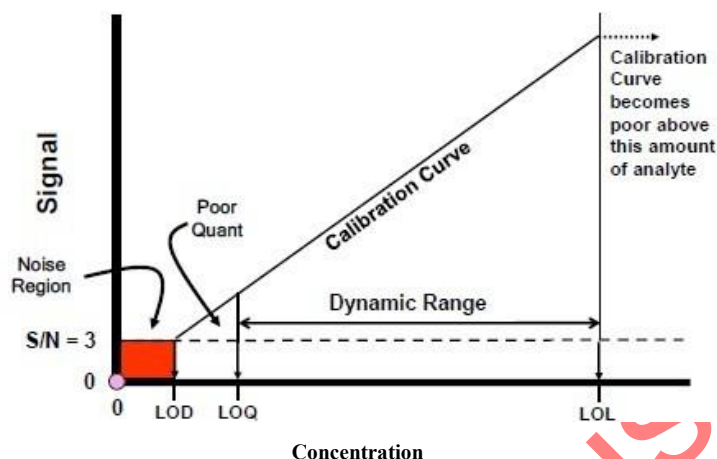


Fig. 3. Calibration curve and figures of merit.

The robustness, or ruggedness, of an analytical method refers to its ability to remain unaffected by small, deliberate variations in experimental conditions, indicating the method's reliability under normal operating conditions [15]. A ruggedness test is a part of method validation and can be considered part of the precision evaluation. Ruggedness is related to repeatability and reproducibility [16]. Robustness is the method's ability to resist change without adapting its initial stable configuration. Ruggedness is a measure of the reproducibility of test results under variations in conditions other than the laboratory.

III. Results and discussion

A. Selectivity and Sensitivity

The selection of wavelengths was based on the normalized intensities for each element achieved for this type of matrix which is cement. This objective was established by running several replicates of the reference materials at the same and varied conditions by standard radial mode and assessing which elements need to be introduced axially (Table 6). Some elements were evaluated using more than one reference wavelength at which energy emission intensities are high and optimum. These are the results when the excited electron returns to the ground state, it can take multiple pathways, each corresponding to the emission of a photon of a specific energy. Thus, one absorption can lead to numerous emissions, therefore, more emission lines at different wavelengths are observed in the atomic spectra. These multiwavelength elements are the Ca (183.801, 220.861, and 616.217), Fe (235.489, 239.562), K (766.490, 769.896), and Na (330.232, 588.995, 589.592). Picking the best wavelength of emission, the actual calibration curves were assessed. Parameters

included high and normalized intensities, goodness of fit (R^2), and standard intensities having a proportional relationship with corresponding concentrations.

All the elements (Ca, Mg, Si, Al, Fe, K, Na, and Ti) show comparable responses except for Na. Na reflects low intensities (denoted with the letter "L") which appreciably vary from one to another. This suggests that Na would give optimum emission intensities if 100% of the sample is nebulized through the plasma via axial mode. Table 7 summarizes the corresponding selective wavelengths and sensitive plasma viewing orientation for analyzing each element of interest.

Table 7. Specific Wavelength and Introduction Mode for Each Element

Element	Ca	Mg	Al	Si	Fe	K	Ti	Na
λ , nm	183.801	280.27	396.153	288.158	239.562	769.896	336.121	330.232
Mode	radial	radial	radial	radial	radial	radial	radial	axial

B. Accuracy and Precision

The oxide constituents of the cement samples were computed (as shown in Table 4). Five replicate runs were made at different conditions (Table 8). The average values were calculated and were presented as the mean \pm standard deviation of these trials. Plotting these results would show the comparison between average percentage recoveries from five replicate runs at different conditions and the target concentrations. These graphs display the profile of the cement types in terms of their stoichiometric oxides. For the Portland type,

Portland Cement Profile

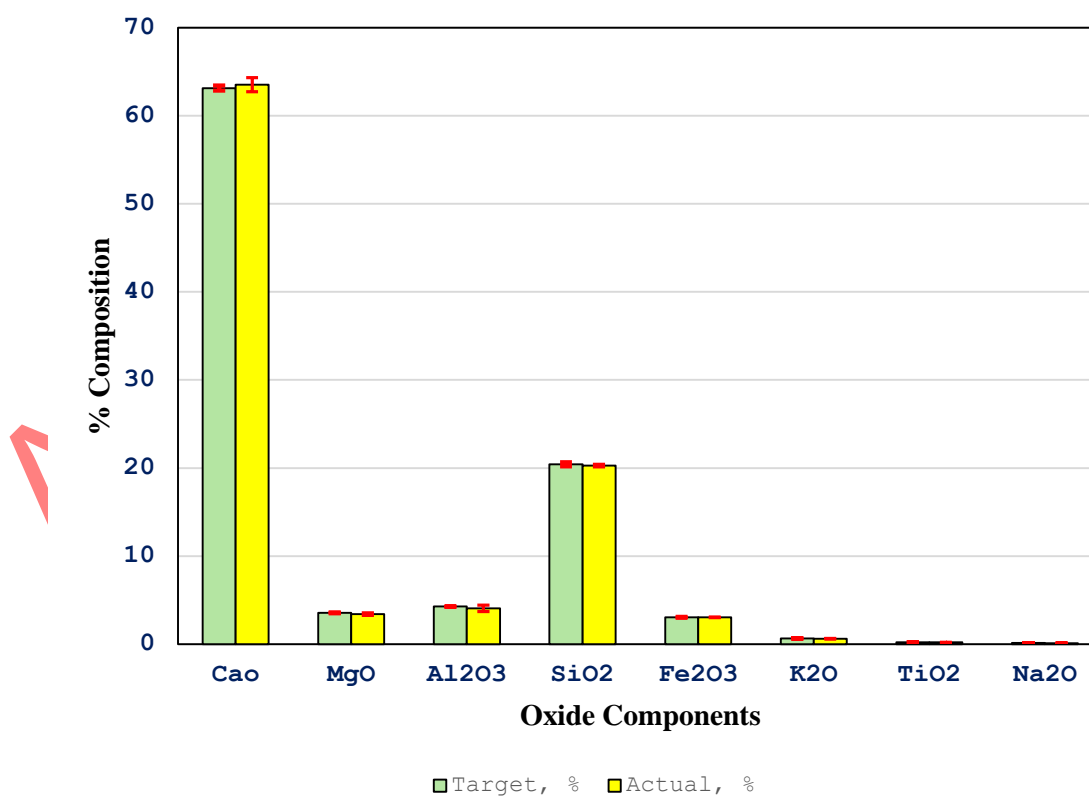


Fig. 4. Oxide constituents' profile of the Portland (SRM 1244b) cement.

For the Blended type (SRM 1881a), the illustration below shows the cement profile and the comparable target and actual percentages with pre-defined and calculated standard deviations, respectively.

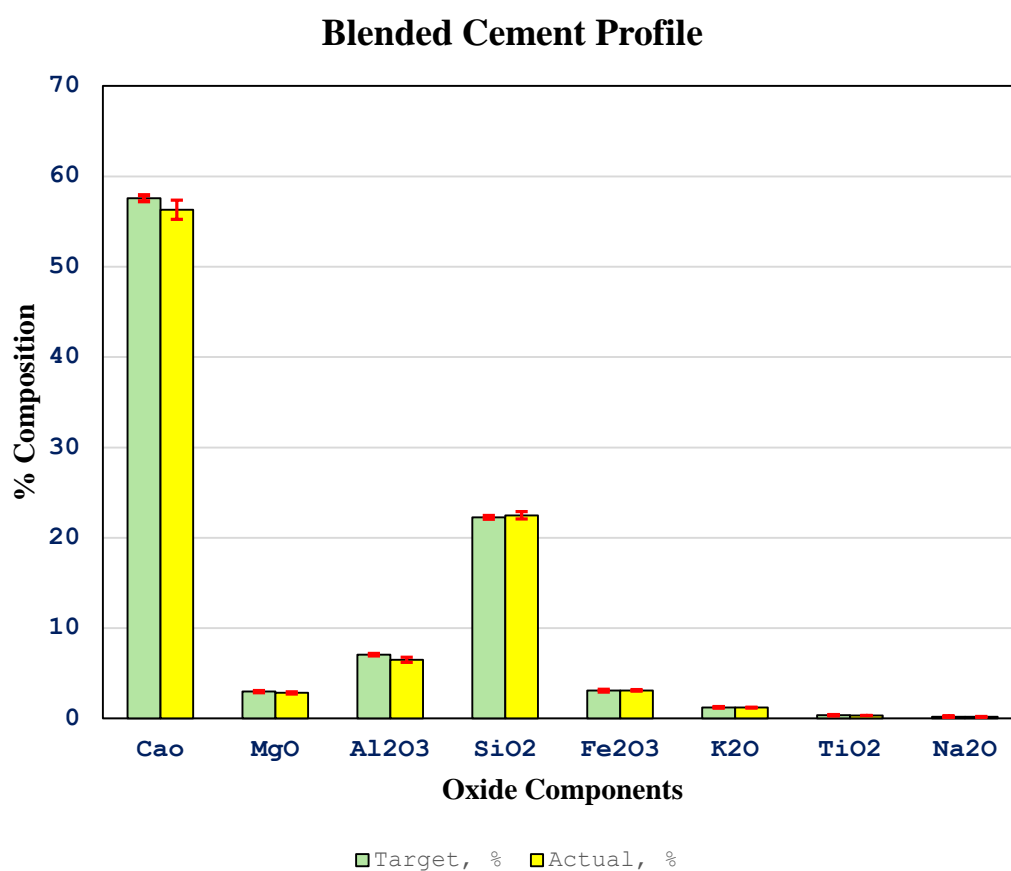


Fig. 9. Oxide constituents' profile of the Blended (SRM 1881a) cement

All components, except Na₂O, have achieved acceptable recovery values at the radial mode of sample introduction. In the case of Na₂O, the axial mode of sample introduction also gave satisfactory recovery results against the targets acclaiming the high accuracy and high precision performance of the closed-system microwave acid digestion and ICP-OES detection.

It can be observed that the corresponding standard deviation values of each component from both cement types were higher than the standard values. This implies that all the uncertainties attributable to the developed method must be established. The estimation utilizes calibration certificates of all paraphernalia involved in the analysis and the factors affecting uncertainty. With this, every effort is made to optimize the measurement procedure. However, the measurement result can only serve as an estimate of the true value, which is inherently unknown. As a result, the closeness of the measured value to the true value cannot be precisely determined, and every measurement carries some degree of uncertainty. The quality and accuracy of a measurement are described by its measurement uncertainty, which defines a range around the measured value within which the true value is expected to lie with a certain probability.

C. Reproducibility Precision

The set of five reproducibility data was evaluated through their HorRat_R values. Concentrations were expressed as a decimal fraction. For reporting HorRat_R values, data were reported as a mass fraction where the units of the numerator and denominator are the same (Table 9).


Table 9. Reproducibility Precision Factors

Constituents	Mass fraction		RSD _R		PRSD _R		HorRat _R	
			Portland	Blend	Portland	Blend	Portland	Blend
CaO	0.6353	0.5631	2.1408	2.1800	2.1429	2.1727	0.9991	1.0034
MgO	0.0342	0.0283	3.3189	3.4138	3.2982	3.3875	1.0063	1.0078
Al ₂ O ₃	0.0406	0.0649	3.2344	3.0141	3.2089	2.9817	1.0079	1.0109
SiO ₂	0.2028	0.2249	2.5407	2.5018	2.5382	2.5055	1.0010	0.9985
Fe ₂ O ₃	0.0307	0.0311	3.3732	3.3662	3.3739	3.3693	0.9998	0.9991
K ₂ O	0.0062	0.0120	4.2831	3.8808	4.2491	3.8695	1.0080	1.0029
TiO ₂	0.0020	0.0033	5.0625	4.7206	4.9945	4.6451	1.0136	1.0163
Na ₂ O	0.0012	0.0017	5.4810	5.1869	5.4068	5.1113	1.0137	1.0148

As expected, the predicted relative standard deviations, PRSD_R, are within 2 to 6 based on the criteria (Table 5). This also holds for the found RSD_R. Since the correspondence between actual and expected is almost exact, the actual HorRat_R values are likely close to 1.0, which means that the method complies with AOAC criteria at high precision, and method reproducibility is acceptable.

D. Linearity

Standard calibration curves and linearity (goodness of fit) were evaluated based on R² (Figure 5).

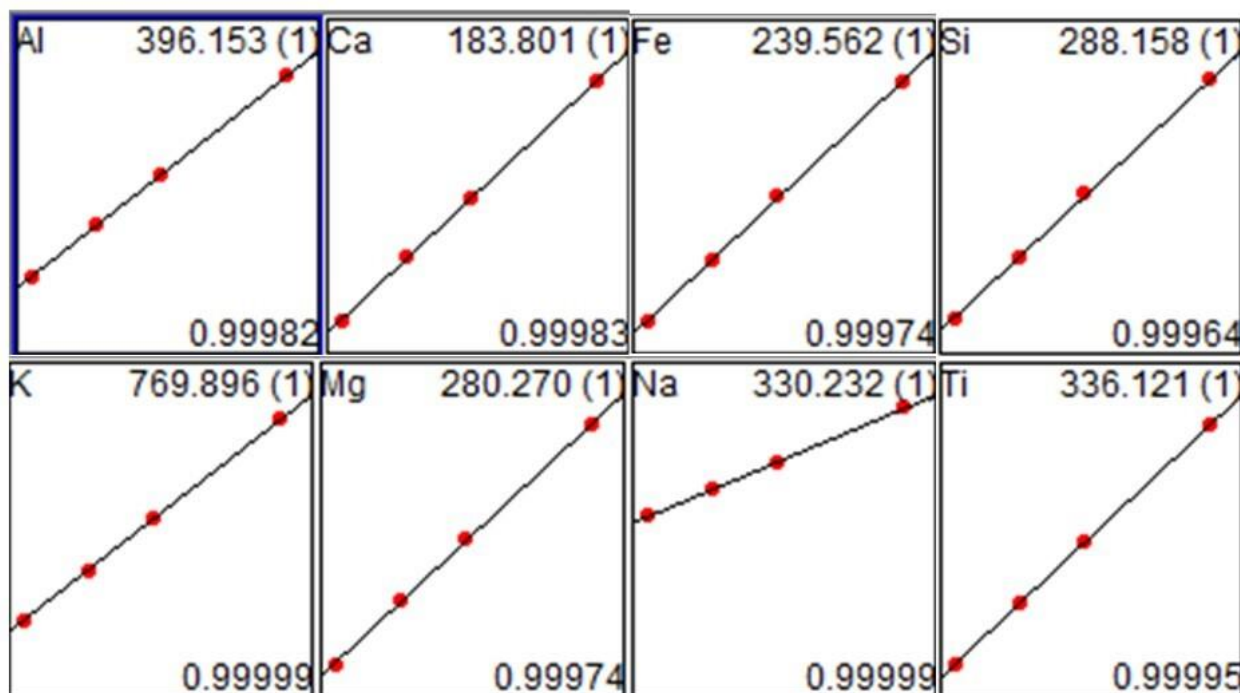

Note: These curves were extracted from the software during the validation run.

Fig. 6. Calibration curves using emission wavelengths of the element of interest.

The R² shows how well a regression model referred to as the series of standard concentrations for each element predicts the outcome of observed data, which is the element of interest in the sample.

In Shimadzu ICP software for ICP-OES, the general form of the calibration equation for a curve fit (often a linear or polynomial regression) in the software is:

$$C = aI^3 + bI^2 + cI + d$$

where, **C** is the calculated concentration of the analyte in the sample.

I is the measured signal intensity.

c is the **slope (sensitivity)** of the linear portion of the calibration curve. It represents the change in concentration per unit of intensity.

d is the **y-intercept (offset)** of the calibration curve. Ideally, this value should be close to zero if the blank measurement is zero and forced through the origin, but it can deviate due to background signals or instrumental drift.

In a simple linear calibration (which is typical for most ICP analyses over a specific range), the **a** and **b** coefficients are typically set to zero, and the equation simplifies to the linear equation,

$$C = cI + d$$

where the values for **c** and **d** are coefficients in the equation used to calculate the concentration from the measured signal intensity [17]. The linearity of the calibration curves was considered suitable for quantitation, with each having a correlation coefficient $R > 0.999$, which gave the correlation of determination, $R^2 > 0.999$ (Table 10).

Table 10. Slope-Intercept Form of the Equation of a Straight Line

Element	WL	R	Order	Origin	<i>c</i>	<i>d</i>	R ²
Al	396.153	0.99982	1	None	0.0049289	-2.912801	0.99964
Ca	183.801	0.99988	1	None	0.2096049	-5.524775	0.99976
Fe	239.562	0.99974	1	None	0.1009963	-6.142189	0.99948
K	769.896	0.99974	1	None	0.0176263	-29.38522	0.99948
Mg	280.270	0.99982	1	None	4.654416e-04	-1.036646	0.99964
Na	330.232	0.99999	1	None	0.0401324	-29.76960	0.99998
Si	288.158	0.99964	1	None	0.0072147	-3.187411	0.99928
Ti	336.121	0.99995	1	None	2.568835e-04	-0.2780649	0.99990

Note: These values were extracted from the software during the validation study.

As observed, not all the lines pass through the origin. This means all the standard intensities have a non-proportional linear relationship with the standard concentrations. Furthermore, the medium used in standards preparation has inherent trace amounts of these elements, resulting in an intrinsic response or matrix effects when running a blank.

E. Range

Samples with high analyte concentrations were diluted with suitable solvents to bring them within the prepared calibration range. The validation has established a dilution pattern (Table 2) to suit the concentration of the analytes present in SRM 1888b and SRM 1881a and fall within the calibration range (Table 3). The evaluation of the method's performance characteristics, particularly the accuracy, precision, and linearity has validated the calibration patterns to be applicable for detecting and quantifying the elemental Al, Ca, Fe, K, Mg, Na, Si, and Ti, and their

equivalent oxide forms when analyzing actual Portland and Blended commercial types of cement. Sample introduction configuration of the plasma was considered a critical part of this method validation study to establish the optimum recoveries of these elements. The dilution step was crucial knowing the analysis would be made simultaneously. However, the linear dynamic range of the ICP when the range of concentrations and their corresponding signals are directly proportional to the concentration of the analyte in the sample was not determined. This can be brought about by preparing a higher standard concentration of each analyte and assessing the goodness of fit of the line to a point when R^2 starts to decrease given its corresponding limit of detection (LOD) and limit of quantification (LOQ).

F. Limit of Detection (LOD)/Limit of Quantification (LOQ)

The limits of detection (LOD) and quantification (LOQ) were determined based on the linearity of the calibration curve and its standard error. LOD is defined as the lowest concentration of an analyte that can be detected by the instrument, and it is not significant to undergo the linearity and precision test. The LOQ is the lowest concentration at which the analyte can be reliably quantified at a 95% confidence interval. Given the calibration curve for each element, the x and y values, and the measured intensities, the residual standard deviation (s_r) for each curve was calculated using the STEYX function in Excel.

Table 11 shows the calculated residual standard deviation for each element with the corresponding LOD and LOQ.

Table 11. LOD and LOQ Values

Elements	Residual Standard Deviation, s_r	Estimated limits (ppm)	
		LOD	LOQ
Ca	0.064	1.0	3.0
Mg	1.5×10^{-5}	0.1	0.3
Si	0.025	12.0	40.0
Al	0.0007	0.5	1.5
K	0.130	25.0	74.0
Fe	0.012	0.4	1.2
Na	0.05	4.1	12.5
Ti	5.4×10^{-5}	0.07	0.21

Sample calculations are shown as part of the supplemental information of this manuscript.

G. Ruggedness/Robustness

The ruggedness of the method was based on the evaluation of the reproducibility precision results. At varying conditions, it was observed that the predicted relative standard deviations ($PRSD_R$) and the relative standard deviations (RSD_R) are within the criteria, resulting in actual $HorRat_R$ values approximately equal to 1.0, meaning the method is rugged.

The method's robustness was challenged by altering the sample size for digestion from 0.2 to 0.5 grams, while maintaining the volumes of the added reagents and detection settings. The objective

was the greater the sample size, the more representative the sample is, and the higher the recovery. However, after the microwave digestion of 0.5 grams samples, the oxides of Ca, Mg, Al, Ti, and Na did not exhibit satisfactory results with increasing sample size. It was noticeable that at 0.5 gram-sample, Na recoveries were no longer good using the axial orientation of the torch. Therefore, the method has decreased its robustness when analyzing greater sample sizes. This could be probably due to salt deposits or eventual clogging inside the nebulizer [18,19]. Digesting 0.5 grams of samples and bringing them up to a certain volume before their introduction into the ICP-OES has resulted in the intensely colored digestates. Potential concentration gradient within the sample solutions could affect the ionization process and varying flow rates or residence times of these solutions within the nebulizer. The likely presence of matrix effects influencing the sensitivity of the analysis would therefore call for additional procedures like standard addition, dilution, separation, filtration, etc., and some of these can greatly contribute to errors in the detection [20,21]. These extra steps can integrate additional uncertainty and variations in the analysis resulting in a bigger bias. Table 12 shows the comparison in terms of test for method robustness between using 0.2- and 0.5-gram samples for analyzing CRMs.

Table 12. Test for Robustness Data

Components	Mode	Targets, %		Actual, %				Recovery, %				NIST compliance			
		Portland	Blended	Portland		Blended		Portland		Blended		Portland		Blended	
				0.2 g	0.5 g	0.2 g	0.5 g	0.2 g	0.5 g	0.2 g	0.5 g	0.2 g	0.5 g	0.2 g	0.5 g
CaO	RADIAL	63.13 ± 0.29	57.58 ± 0.34	63.1	55.83	57.58	57.79	99.96	88.4	100.0	100.4	✓		✓	✓
MgO		3.562 ± 0.057	2.981 ± 0.077	3.3	----	2.8	2.7	93.6	----	94.6	90.7	✓		✓	✓
Al ₂ O ₃		4.277 ± 0.036	7.06 ± 0.081	4.1	4.0	6.5	4.8	95.9	93.2	92.1	68.9	✓	✓	✓	
SiO ₂		20.42 ± 0.23	22.26 ± 0.15	20.5	19.1	22.9	20.4	100.3	93.7	102.8	91.5	✓	✓	✓	✓
Fe ₂ O ₃		3.062 ± 0.053	3.09 ± 0.11	3.06	2.82	3.23	3.13	99.9	92.1	104.5	101.3	✓	✓	✓	✓
K ₂ O		0.658 ± 0.017	1.228 ± 0.029	0.64	0.63	1.23	1.18	97.3	95.7	100	96.03	✓	✓	✓	✓
TiO ₂		0.2316 ± 0.0076	0.3663 ± 0.0030	0.23	0.19	0.33	0.30	98.7	82.9	91.1	82.9	✓		✓	
Na ₂ O	AXIAL	0.1364 ± 0.0044	0.199 ± 0.007	0.13	0.11	0.18	0.17	93.8	80.7	91.0	85.4	✓		✓	

Note: The analysis of magnesium using 0.5-gram samples has exceeded the equipment detection limit. Potential saturation during detection has occurred, signifying response optimization through dilution or decrease in sample size. Other metals like Ca, Al, Ti, and Na have reached the saturation level of detection and thus affected the linear range, accordingly. Consequently, this has resulted in slightly lowered intensities that caused slightly lowered concentrations. The (✓) indicates assay conformity of each component according to the NIST Certificates of Analysis for Portland (SRM 1888b) and Blended (SRM 1881a).

Conclusions and Recommendations

The findings of this study demonstrate that the developed analytical method is reliable and effective for both Portland and blended cement samples. With appropriate in-house validation using local reference materials, the protocol is also expected to apply to masonry-type samples. The method's performance was supported by acceptable figures of merit obtained from 0.5-gram samples digested through a closed-vessel microwave system and analyzed following the established calibration pattern. Routine operations using microwave-ICP were successfully carried out with dedicated polypropylene or non-contaminating laboratory ware to minimize sodium, silicon, and aluminum carryover. Optimal plasma viewing modes were identified, with axial mode required for sodium and radial torch alignment for aluminum, calcium, magnesium, potassium, silicon, titanium, and iron determinations.

Despite the overall success of the validation study, the method showed limitations when applied to larger (1.0-gram) samples, indicating a lack of robustness at higher sample masses. Future work should include establishing the instrument's linear dynamic range by verifying calibration performance (e.g., R^2 values) across extended concentration levels. Routine quality control measures such as periodic introduction of check samples or distilled water, regular inspection of the sample introduction system, and proper maintenance of pump tubing and rinse solutions are essential to ensuring analytical integrity. Finally, estimating the measurement uncertainty is recommended to better quantify potential sources of error and strengthen confidence in the method's applicability for routine cement analysis.

Conflicts of Interest

The author declares that no conflicts of interest exist.

List of Nomenclature

C	concentration of the analyte
HorRat_R	Horwitz ratio (HorRat) or a normalized performance parameter indicating the acceptability of the method of analysis with respect to reproducibility precision.
PRSD_R	predicted reproducibility relative standard deviation (%)
RSD_R	reproducibility relative standard deviation (%)
s_f	Residual Standard Deviation

References

1. José N, Ahmed H, Miguel B, Luís E, Jorge de B. Magnesia (MgO) Production and Characterization, and Its Influence on the Performance of Cementitious Materials: A Review. *Materials* [Internet]. 2020 Oct 23;13(21):4752. Available from:

<https://www.ncbi.nlm.nih.gov/pmc/articles/PMC7660681/>.

<https://doi.org/10.3390/ma13214752>

٢. ASTM C114-00. Standard Test Methods for Chemical Analysis of Hydraulic Cement.
٣. Bediako M, Amankwah EO. Analysis of Chemical Composition of Portland Cement in Ghana: A Key to Understand the Behavior of Cement. *Advances in Materials Science and Engineering*. 2015;2015:1–5. <https://doi.org/10.1155/2015/349401>
٤. Gustavo González A, Ángeles Herrador M. A practical guide to analytical method validation, including measurement uncertainty and accuracy profiles. *TrAC Trends in Analytical Chemistry* [Internet]. 2007 Mar;26(3):227–38. Available from: <https://www.sciencedirect.com/science/article/pii/S0165993607000118#>.
<https://doi.org/10.1016/j.trac.2007.01.009>
٥. Nayak A, Parui K, Sharma S, Ratha S. Study and Analysis of Atomic Spectra. *International Journal of Scientific and Research Publications (IJSRP)*. 2020 Nov 24;10(11):946–55. <https://doi.org/10.29322/ijsrp.10.11.2020.p10787>
٦. Hou X, Amais RS, Jones BT, Donati GL. Inductively Coupled Plasma Optical Emission Spectrometry. *Encyclopedia of Analytical Chemistry*. 2016 Mar 9;1–25. <https://doi.org/10.1002/9780470027318.a5110.pub3>
٧. Sullivan MP, Holtkamp HU, Meier SM, Hartinger CG. The Analysis of Therapeutic Metal Complexes and Their Biomolecular Interactions. Elsevier eBooks. 2017 Jan 1;355–86. <https://doi.org/10.1016/b978-0-12-803814-7.00010-1>
٨. Iti S. ICP-OES: An Advance Tool in Biological Research. *Open Journal of Environmental Biology*. 2020 Dec 26;027–33. <https://doi.org/10.17352/ojeb.000018>
٩. Betz JM, Brown PN, Roman MC. Accuracy, precision, and reliability of chemical measurements in natural products research. *Fitoterapia* [Internet]. 2011 Jan;82(1):44–52. Available from: <https://www.ncbi.nlm.nih.gov/pmc/articles/PMC3026088/>.
<https://doi.org/10.1016/j.fitote.2010.09.011>
١٠. Horwitz W. The Variability of AOAC Methods of Analysis as Used in Analytical Pharmaceutical Chemistry. *Journal of AOAC INTERNATIONAL*. 1977 Nov 1;60(6):1355–63. <https://doi.org/10.1093/jaoac/60.6.1355>
١١. McClure FD, Lee JK. Computation of HORRAT Values. *Journal of AOAC International*. 2003 Sep 1;86(5):1056–8. <https://doi.org/10.1093/jaoac/86.5.1056>
١٢. AOAC Guidelines for Single-Laboratory Validation of Chemical Methods for Dietary Supplements and Botanicals (2012) *Official Methods of Analysis*, 19th Ed., Appendix K, AOAC INTERNATIONAL, Gaithersburg, MD.
١٣. Moosavi SM, Ghassabian S. Linearity of Calibration Curves for Analytical Methods: A Review of Criteria for Assessment of Method Reliability. *Calibration and Validation of Analytical Methods - A Sampling of Current Approaches* [Internet]. 2018 Apr 25; Available from: <https://www.intechopen.com/books/calibration-and-validation-of->

[analytical-methods-a-sampling-of-current-approaches/linearity-of-calibration-curves-for-analytical-methods-a-review-of-criteria-for-assessment-of-method.](https://doi.org/10.5772/intechopen.72932)

<https://doi.org/10.5772/intechopen.72932>

14. Al-Hakkani MF. A rapid, developed and validated RP-HPLC method for determination of azithromycin. *SN Applied Sciences*. 2019 Feb 11;1(3).
15. Tome T, Žigart N, Časar Z, Obreza A. Development and Optimization of Liquid Chromatography Analytical Methods by Using AQbD Principles: Overview and Recent Advances. *Organic Process Research & Development*. 2019 Aug 14;23(9):1784–802. <https://doi.org/10.1021/acs.oprd.9b00238>
16. Vander Heyden Y, Massart DL. Chapter 3 Review of the use of robustness and ruggedness in analytical chemistry. *Data Handling in Science and Technology*. 1996;79–147. [https://doi.org/10.1016/s0922-3487\(96\)80016-5](https://doi.org/10.1016/s0922-3487(96)80016-5)
17. Chen HY, Chen C. Evaluation of Calibration Equations by Using Regression Analysis: An Example of Chemical Analysis. *Sensors* [Internet]. 2022 Jan 7;22(2):447. Available from: <https://pmc.ncbi.nlm.nih.gov/articles/PMC8780995/>
18. Clarice D.B. Amaral, Fialho LL, Felipe P.R. Camargo, Camillo Pirola, Nóbrega JA. Investigation of analyte losses using microwave-assisted sample digestion and closed vessels with venting. *Talanta*. 2016 Nov 1;160:354–9.
19. Genovesi L. *Considerations for optimizing an ICP sample introduction system*. 2017 Dec 15.
20. Barros AI, Pinheiro FC, Nóbrega JA. Calibration strategies to correct for matrix effects in direct analysis of urine by ICP-OES: internal standardization and multi-energy calibration. *Analytical Methods*. 2019;11(27):3401-3409. doi:10.1039/C9AY00907H.
21. Karlsson S, Andersson JT, Petersson P. Matrix effects in simultaneous microwave-induced plasma optical emission spectrometry: new perspectives on an old problem. *Journal of Analytical Atomic Spectrometry*. 2023;38(9):1983-1992. doi:10.1039/D3JA00061C.

Supplemental Information

Sample calculations for the method's LOD and LOQ are shown below. For example, in the case of Na, using the STEYX function of EXCEL, the residual standard deviation of the calibration curve was calculated as $s_r = 0.05$. Given the slope, $S = 0.0401324$, the estimated LOD was found to be:

$$LOD = 3.3 \times s_r / S$$

$$LOD = 3.3 \times \frac{0.05}{0.0401324}$$

$$LOD \cong 4.1 \text{ ppm}$$

Thus, the LOQ:

$$LOQ = 10 \times s_r / S$$

$$LOQ = 10 \times \frac{0.05}{0.0401324}$$

$$LOQ \cong 12.5 \text{ ppm}$$

Another example is for Mg, given $s_r = 1.5 \times 10^{-5}$ with a slope $S = 4.6544 \times 10^{-4}$,

$$LOD = 3.3 \times \frac{0.000015}{0.00046544}$$

$$LOD \cong 0.1 \text{ ppm}$$

Thus, the LOQ:

$$LOQ = 10 \times s_r / S$$

$$LOQ = 10 \times \frac{0.000015}{0.00046544}$$

$$LOQ \cong 0.3 \text{ ppm}$$

Table 6. Recorded Intensities for Ca, Mg, Si, Al, Fe, K, Na and Ti

Sample Name	Al 396.153	Ca 183.801	Ca 220.861	Ca 616.217	Fe 235.489	Fe 239.562	K 766.490	K 769.896	Mg 280.270	Na 330.232	Na 588.995	Na 589.592	Si 288.158	Ti 336.121
Portland 1	12100	371000	323000	414000	26500	22100	5310	5540	!!!!!!!	2760 L	181	231	107000	1270
Portland 2	14300	375000	323000	418000	26700	21100	5310	5520	!!!!!!!	2780 L	184	234	112000	1260
Portland 3	14700	378000	321000	413000	26500	21300	5270	5530	!!!!!!!	2840 L	187	235	106000	1250
Portland 4	14800	377000	323000	414000	26400	21200	5250	5530	!!!!!!!	2820 L	179	227	109000	1250
Portland 5	14700	376000	324000	411000	26400	21500	5270	5540	!!!!!!!	2840 L	182	231	109000	1250
Portland 6	15600	379000	323000	420000	26600	21200	5290	5540	!!!!!!!	2840 L	181	230	108000	1260
Portland 7	15300	377000	321000	415000	26500	21000	5270	5520	!!!!!!!	2830 L	180	229	112000	1250
Mean Average	14500	376143	322571	415000	26514	21343	5281	5531	#DIV/0!	#DIV/0!	182	231	109000	1256
Standard Dev.	1142	2610	1134	3055	107	369	23	9	#DIV/0!	#DIV/0!	3	3	2309	8
%RSD	7.873	0.694	0.352	0.736	0.403	1.729	0.429	0.163	#DIV/0!	#DIV/0!	1.488	1.199	2.119	0.627
Portland 8	16100	337000	294000	364000	25800	21900	10300	10600	!!!!!!!	2850 L	251	303	119000	1990
Portland 9	20200	346000	291000	373000	25800	21100	10300	10500	!!!!!!!	2840 L	254	305	121000	2000
Portland 10	17900	337000	291000	367000	25700	21500	10300	10500	!!!!!!!	2840 L	253	306	118000	1990
Portland 11	19500	335000	293000	367000	25600	21400	10300	10500	!!!!!!!	2830 L	254	307	118000	1980
Portland 12	23300	345000	297000	374000	25700	21500	10300	10600	!!!!!!!	2830 L	250	303	116000	1990
Portland 13	20000	348000	290000	376000	25900	21200	10300	10500	!!!!!!!	2810 L	258	309	120000	2010
Portland 14	21000	343000	292000	371000	25600	21500	10200	10500	!!!!!!!	2820 L	249	301	112000	1990
Mean Average	19714	341571	292571	370286	25729	21443	10286	10529	#DIV/0!	#DIV/0!	253	305	117714	1993
Standard Dev.	2281	5159	2370	4386	111	257	38	49	#DIV/0!	#DIV/0!	3	3	2984	10
%RSD	11.572	1.510	0.810	1.185	0.432	1.200	0.367	0.463	#DIV/0!	#DIV/0!	1.203	0.897	2.535	0.477
Portland 15	15400 L	467000	403000	447000	32100	24800	5290	4450	22100	6300 L	271	521	137000	1240
Portland 16	18000	473000	401000	454000	32100	24400	5220	4280	22200	6190 L	245	497	124000	1240
Portland 17	18600	468000	399000	450000	32000	24700	4540	3470	22100	5920 L	166	422	115000	1230
Portland 18	18600	471000	405000	449000	31800	24700	5130	4040	22100	6050 L	311	569	113000	1220
Portland 19	18600	474000	399000	457000	32100	24700	5210	3780	22300	5650 L	293	563	115000	1240
Portland 20	19500	479000	400000	457000	32000	24500	5200	3820	22200	5590 L	262	530	109000	1240
Portland 21	19100	471000	399000	451000	31700	24700	4660	3380	22000	5520 L	159	425	107000	1230
Mean Average	18257	471857	400857	452143	31971	24643	5036	3889	22143	#DIV/0!	244	504	117143	1234
Standard Dev.	1344	4018	2340	3934	160	140	303	396	98	#DIV/0!	60	60	10303	8

%RSD	7.361	0.851	0.584	0.870	0.502	0.567	6.022	10.188	0.441	#DIV/0!	24.406	11.932	8.795	0.637
Blended 1	19900 L	419000	357000	403000	30300	24500	9500	7910	18400	4610 L	83.7	351	113000	2090
Blended 2	24300	432000	362000	407000	30200	24300	9840	8460	18500	4960 L	188	449	109000	2070
Blended 3	21900	426000	359000	402000	30100	24400	9910	8450	18400	4880 L	205	466	110000	2070
Blended 4	23700	425000	358000	402000	30000	24300	10300	8760	18400	4800 L	302	560	110000	2060
Blended 5	27800	436000	360000	413000	30200	24300	9570	7880	18700	4640 L	119	388	109000	2090
Blended 6	24000	428000	352000	405000	30000	24000	9650	8160	18500	4610 L	156	418	111000	2070
Blended 7	25400	422000	363000	404000	30100	24400	10000	8570	18400	4710 L	213	472	105000	2060
Mean Average	23857	426857	358714	405143	30129	24314	9824	8313	18471	#DIV/0!	181	443	109571	2073
Standard Dev.	2505	5786	3638	3891	111	157	278	337	111	#DIV/0!	71	67	2440	13
%RSD	10.501	1.355	1.014	0.960	0.369	0.647	2.834	4.048	0.602	#DIV/0!	39.223	15.186	2.227	0.605
Blended 8	15000	314000	274000	339000	24000	20400	9640	9880		2650 L	233	283	111000	1860
Blended 9	18600	319000	269000	344000	23800	19500	9500	9720	!!!!!!	2620 L	234	282	111000	1850
Blended 10	16500	310000	268000	339000	23700	19800	9460	9690	!!!!!!	2610 L	233	282	109000	1840
Blended 11	18100	310000	272000	341000	23800	19800	9510	9730	!!!!!!	2630 L	236	285	110000	1840
Blended 12	21700	322000	277000	349000	24000	20000	9620	9860	!!!!!!	2640 L	233	282	108000	1860
Blended 13	18500	323000	269000	348000	24000	19700	9550	9760	!!!!!!	2610 L	239	287	111000	1870
Blended 14	19600	319000	272000	345000	23900	20000	9510	9760	!!!!!!	2630 L	232	280	105000	1850
Mean Average	18286	316714	271571	343571	23886	19886	9541	9771	#DIV/0!	#DIV/0!	234	283	109286	1853
Standard Dev.	2143	5407	3207	4077	121	285	66	72	#DIV/0!	#DIV/0!	2	2	2215	11
%RSD	11.718	1.707	1.181	1.187	0.509	1.435	0.694	0.734	#DIV/0!	#DIV/0!	1.037	0.816	2.026	0.601
Blended 15	18600 L	391000	333000	376000	28300	22800	8850	7370	17100	4290 L	78	327	105000	1950
Blended 16	22500	398000	334000	376000	27900	22400	9080	7800	17100	4580 L	173	415	101000	1910
Blended 17	20200	392000	331000	371000	27800	22400	9130	7780	16900	4490 L	189	430	101000	1910
Blended 18	22000	395000	332000	373000	27900	22500	9520	8120	17100	4450 L	280	520	102000	1910
Blended 19	26000	406000	336000	385000	28200	22600	8930	7350	17500	4330 L	111	362	101000	1950
Blended 20	22200	397000	326000	376000	27800	22300	8960	7570	17100	4270 L	144	388	103000	1920
Blended 21	23600	393000	338000	377000	28000	22800	9310	7980	17200	4380 L	199	439	98200	1920
Mean Average	22157	396000	332857	376286	27986	22543	9111	7710	17143	#DIV/0!	168	412	101600	1924
Standard Dev.	2362	5099	3848	4386	195	199	235	294	181	#DIV/0!	66	62	2094	18
%RSD	10.661	1.288	1.156	1.166	0.697	0.882	2.581	3.814	1.057	#DIV/0!	39.171	15.052	2.061	0.942

Note: Samples were aspirated using the default sample introduction mode which was the radial orientation of the nebulizer to the torch being the plasma source.

Table 8. Summary of Concentrations of Oxide Constituents Replicated at Different Conditions and Corresponding Percentage Recovery (R)

Oxide Component s	Mode	Target, %		Actual, %										Average, %		Recovery (R), %	
		Portlan d	Blende d	Portland					Blended					Portlan d	Blende d	Portlan d	Blende d
				1 st	2 nd	3 rd	4 th	5 th	1 st	2 nd	3 rd	4 th	5 th				
CaO	Radial	63.13 ± 0.29	57.58 ± 0.34	62.64	64.22	63.02	64.52	63.26	57.19	57.23	55.55	54.87	56.69	63.53 ± 0.80	56.31 ± 1.05	100.6	97.8
MgO		3.562 ± 0.057	2.981 ± 0.077	3.27	3.38	3.31	3.45	3.67	2.79	2.78	2.73	2.79	3.06	3.42 ± 0.16	2.83 ± 0.13	95.9	95.0
Al ₂ O ₃		4.277 ± 0.036	7.06 ± 0.081	4.02	3.52	3.96	4.30	4.48	6.23	6.19	6.77	6.56	6.73	4.06 ± 0.37	6.50 ± 0.27	94.9	93.1
SiO ₂		20.42 ± 0.23	22.26 ± 0.15	20.10	20.23	20.27	20.30	20.52	22.15	23.16	22.22	22.51	22.38	20.28 ± 0.15	22.48 ± 0.40	99.3	101.0
Fe ₂ O ₃		3.062 ± 0.053	3.09 ± 0.11	3.12	3.07	3.07	3.02	3.05	3.28	3.10	3.07	3.03	3.07	3.07 ± 0.04	3.11 ± 0.10	100.1	100.6
K ₂ O		0.658 ± 0.017	1.228 ± 0.029	0.63	0.64	0.60	0.60	0.67	1.22	1.22	1.15	1.15	1.27	0.63 ± 0.03	1.20 ± 0.05	94.8	98.1
TiO ₂		0.2316 ± 0.0076	0.3663 ± 0.0030	0.20	0.20	0.22	0.19	0.21	0.33	0.33	0.32	0.33	0.33	0.20 ± 0.01	0.33 ± 0.004	91.4	89.8
Na ₂ O	Axial	0.1364 ± 0.0044	0.199 ± 0.007	0.12	0.13	0.10	0.09	0.16	0.19	0.13	0.18	0.13	0.25	0.12 ± 0.03	0.18 ± 0.05	91.3	90.7

Note: The average recoveries are reported as mean ± standard deviation of the five reproducible replicates. Oxides of Ca, Mg, Al, Si, Fe, K, and Ti have achieved acceptable recovery values radially. The axial mode of sample introduction on Na gave satisfactory recoveries for both cement types.